

Day-Bailes tapered locking suction tube for micro and spine surgery

Designed by Dr. John Diaz Day, Associate Professor of Neurosurgery at the University of Texas, San Antonio and Dr. Julian E. Bailes, Chairman and Chief of the Department of Neurosurgery at West Virginia University

Ergonomic and well balanced, with precise intraoperative suction control


Designed for a wide range of delicate spine and microsurgical procedures. This versatile system combines an ergonomic, surgeon-friendly design with the freedom and flexibility of 40 interchangeable tube options. The slim, lightweight handpiece reduces fatigue during long procedures, while providing precise, intraoperative suction control with thumb or forefinger.

Economy and convenience

Interchangeable, quick-release tubes come in four lengths, ten diameters and in one French increments, allowing instrument customization to fit your preferences and procedure. Unlike with other models, tubes may be easily changed out during surgery without discarding the handpiece.

Day-Bailes suction tube advantages

- Precisely milled compression fitting and guide pin align to securely lock and hold the tube in place
- Inner and outer cannula tapered to reduce clogging
- Malleable tubes come pre-bent to 30° to prevent kinking
- Short, medium, long and extra long, in sizes 3 to 12 French
- Your choice of flat or triangular handpiece
- Tubes may be ordered separately to minimize replacement costs

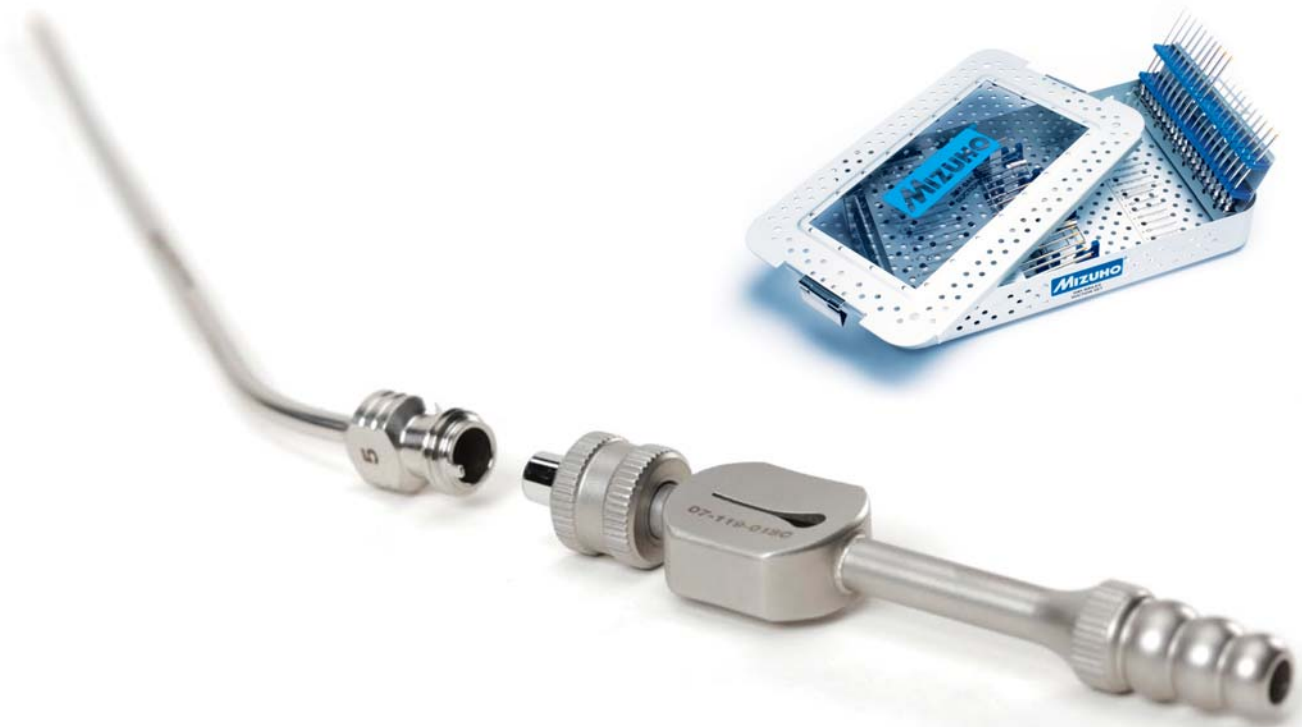



Kelly endonasal suction tubes



Atraumatic suction tubes



Triangular and flat handpieces

Product information

Long and extra long tubes feature a gold distal tip to determine if the tube has been drilled down and for intraoperative orientation.

Tube lengths from distal tip to the bend:

Short: 7.6 cm Medium: 9.5 cm
 Long: 11.4 cm Extra long: 13.3 cm

Flat handpiece option lengths:

Short: 17.5 cm Medium: 19.5 cm
 Long: 21.7 cm Extra long: 23.7 cm

Triangular handpiece option lengths:

Short: 20.5 cm Medium: 22.2 cm
 Long: 24.6 cm Extra long: 26.5 cm

Ordering Information

Catalog No. Description

Handles

07-119-01SC Flat teardrop suction handle w/safety coupling
 07-DAY-01SC Day suction tube handle w/safety coupling

Tubes

07-DAY-03SSC Day-Bailes suction 3 French short locking
 07-DAY-03MSC Day-Bailes suction 3 French medium locking
 07-DAY-03LSC Day-Bailes suction 3 French long locking
 07-DAY-03XLSC Day-Bailes suction 3 French x-long locking
 07-DAY-04SSC Day-Bailes suction 4 French short locking
 07-DAY-04MSC Day-Bailes suction 4 French medium locking
 07-DAY-04LSC Day-Bailes suction 4 French long locking
 07-DAY-04XLSC Day-Bailes suction 4 French x-long locking suction
 07-DAY-05SSC Day-Bailes suction 5 French short locking
 07-DAY-05MSC Day-Bailes suction 5 French medium locking
 07-DAY-05LSC Day-Bailes suction 5 French long locking
 07-DAY-05XLSC Day-Bailes suction 5 French x-long locking suction
 07-DAY-06SSC Day-Bailes suction 6 French short locking
 07-DAY-06MSC Day-Bailes suction 6 French medium locking
 07-DAY-06LSC Day-Bailes suction 6 French long locking
 07-DAY-06XLSC Day-Bailes suction 6 French x-long locking
 07-DAY-07SSC Day-Bailes suction 7 French short locking
 07-DAY-07MSC Day-Bailes suction 7 French medium locking
 07-DAY-07LSC Day-Bailes suction 7 French long locking
 07-DAY-07XLSC Day-Bailes suction 7 French x-long locking
 07-DAY-08SSC Day-Bailes suction 8 French short locking
 07-DAY-08MSC Day-Bailes suction 8 French medium locking

Catalog No. Description

Tubes

07-DAY-08LSC Day-Bailes suction 8 French long locking
 07-DAY-08XLSC Day-Bailes suction 8 French x-long locking
 07-DAY-09SSC Day-Bailes suction 9 French short locking
 07-DAY-09MSC Day-Bailes suction 9 French medium locking
 07-DAY-09LSC Day-Bailes suction 9 French long locking
 07-DAY-09XLSC Day-Bailes suction 9 French x-long locking
 07-DAY-10SSC Day-Bailes suction 10 French short locking
 07-DAY-10MSC Day-Bailes suction 10 French medium locking
 07-DAY-10LSC Day-Bailes suction 10 French long locking
 07-DAY-10XLSC Day-Bailes suction 10 French x-long locking
 07-DAY-11SSC Day-Bailes suction 11 French short locking
 07-DAY-11MSC Day-Bailes suction 11 French medium locking
 07-DAY-11LSC Day-Bailes suction 11 French long locking
 07-DAY-11XLSC Day-Bailes suction 11 French x-long locking
 07-DAY-12SSC Day-Bailes suction 12 French short locking
 07-DAY-12MSC Day-Bailes suction 12 French medium locking
 07-DAY-12LSC Day-Bailes suction 12 French long locking
 07-DAY-12XLSC Day-Bailes suction 12 French x-long locking

Optional Suction Tubes

07-961-05L Multi hole atraumatic suction tube 5 French long
 07-961-05XL Multi hole atraumatic suction tube 5 French XL
 07-961-07L Multi hole atraumatic suction tube 7 French long
 07-961-07XL Multi hole atraumatic suction tube 7 French XL
 07-961-09L Multi hole atraumatic suction tube 9 French long
 07-961-09XL Multi hole atraumatic suction tube 9 French XL
 07-961-12L Multi hole atraumatic suction tube 12 French long
 07-961-12XL Multi hole atraumatic suction tube 12 French XL
 07-962-05D Kelly endonasal suction tube 5 French 150mm Down
 07-962-05U Kelly endonasal suction tube 5 French 150mm Up
 07-962-07D Kelly endonasal suction tube 7 French 150mm Down
 07-962-07U Kelly endonasal suction tube 7 French 150mm Up

Mizuho America, Inc.
 133 Brimbal Avenue
 Beverly, MA 01915
 mizuho.com

Call Toll-free:: 800 669 2547
 email:: mizuho@mizuho.com